

Dr. Zaghi: Light Scalpel Frenectomy Settings



4 WATTS


Repeat 29 Hz F2-2

10 msec
40 mJ
1160 mW Av
29.0 % Duty

MUCOSA MODE

SUPER PULSE

Mucosa

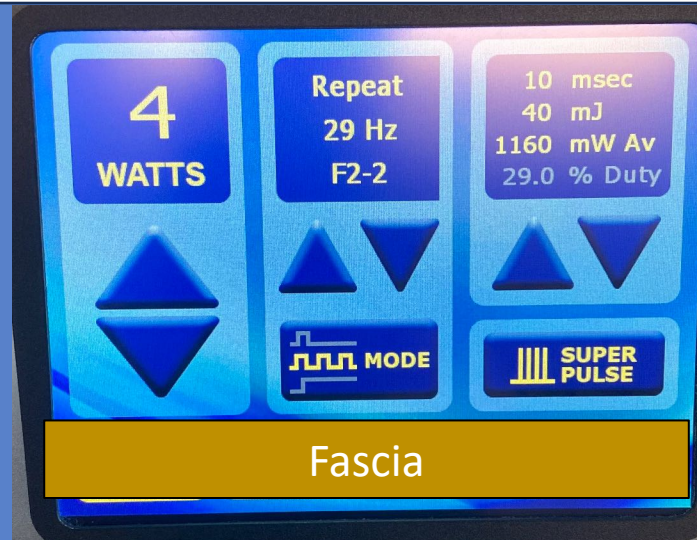

4 WATTS


Repeat 29 Hz F2-2

10 msec
40 mJ
1160 mW Av
29.0 % Duty

MUCOSA MODE

SUPER PULSE

Fascia

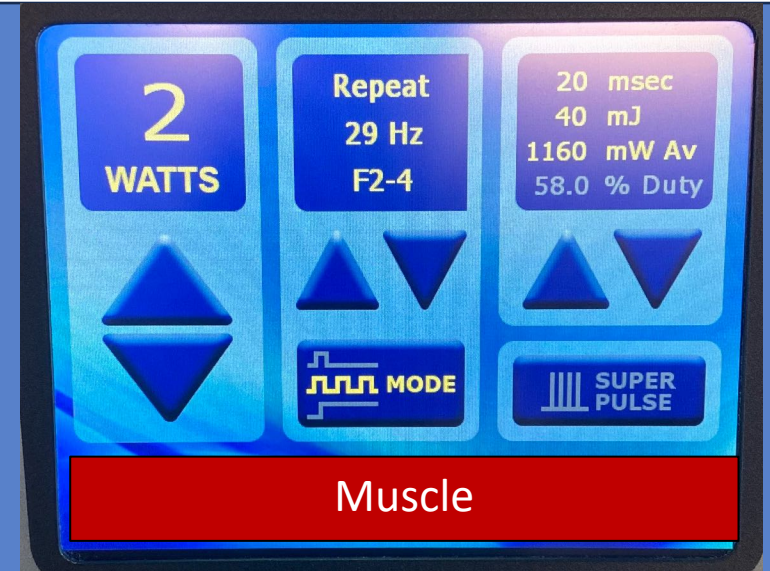

2 WATTS


Repeat 29 Hz F2-4

20 msec
40 mJ
1160 mW Av
58.0 % Duty

MUCOSA MODE

SUPER PULSE

Muscle

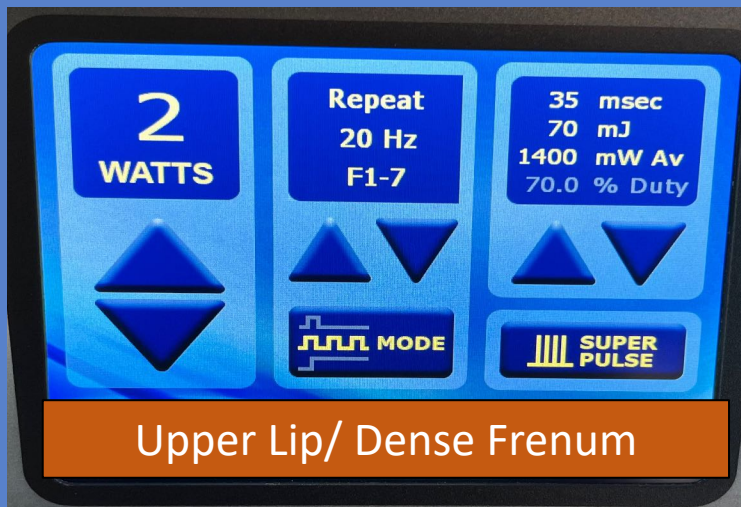

2 WATTS


Repeat 20 Hz F1-7

35 msec
70 mJ
1400 mW Av
70.0 % Duty

MUCOSA MODE

SUPER PULSE

Upper Lip/ Dense Frenum

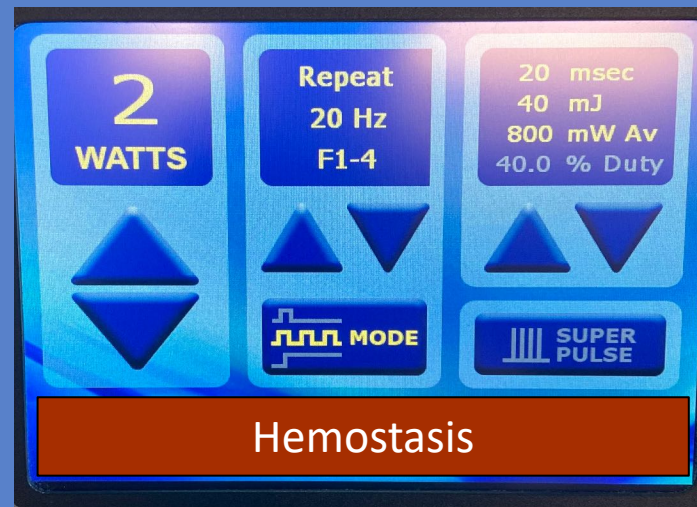

2 WATTS


Repeat 20 Hz F1-4

20 msec
40 mJ
800 mW Av
40.0 % Duty

MUCOSA MODE

SUPER PULSE

Hemostasis

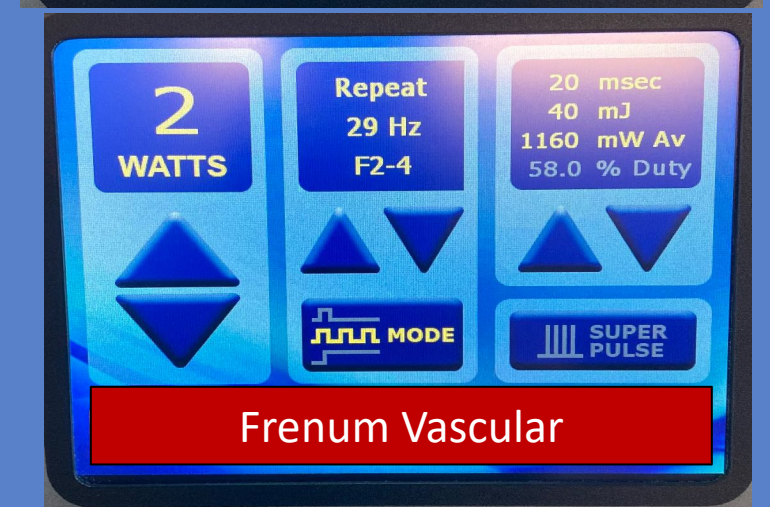

2 WATTS

Repeat 29 Hz F2-4

20 msec
40 mJ
1160 mW Av
58.0 % Duty

MUCOSA MODE


SUPER PULSE

Frenum Vascular

Non-Ablative Tissue Remodeling Settings



Non-Ablative: *Recommended*

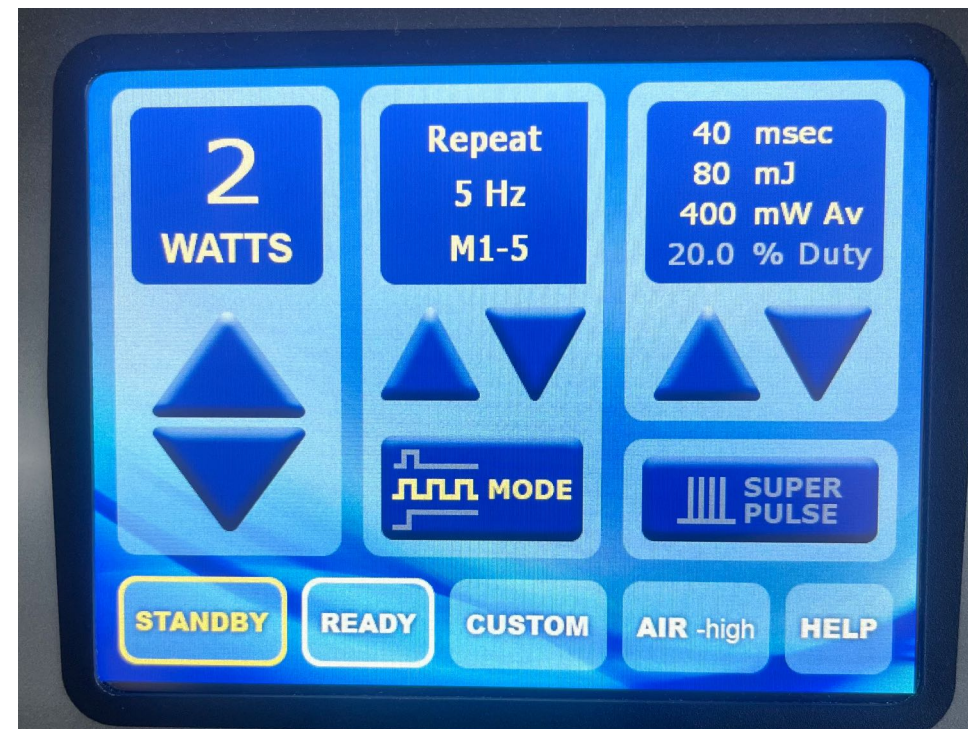

Non-Ablative: Slower and Less Intense